

INSTALLATION INSTRUCTIONS
for the
REPLENISHMENT TIMER
on the
***Kodak X-Omat* M35A-M and M35-M PROCESSORS**
Service Codes: 3133, 3408

***Kodak Professional Industrial* M35I-M PROCESSOR**
Service Code: 3400

***Kodak Professional Industrex* M35 PROCESSOR**
Service Code: 3414



HEALTH IMAGING



PLEASE NOTE The information contained herein is based on the experience and knowledge relating to the subject matter gained by Eastman Kodak Company prior to publication.

No patent license is granted by this information.

Eastman Kodak Company reserves the right to change this information without notice, and makes no warranty, express or implied, with respect to this information. Kodak shall not be liable for any loss or damage, including consequential or special damages, resulting from any use of this information, even if loss or damage is caused by Kodak's negligence or other fault.



This equipment includes parts and assemblies sensitive to damage from electrostatic discharge. Use caution to prevent damage during all service procedures.

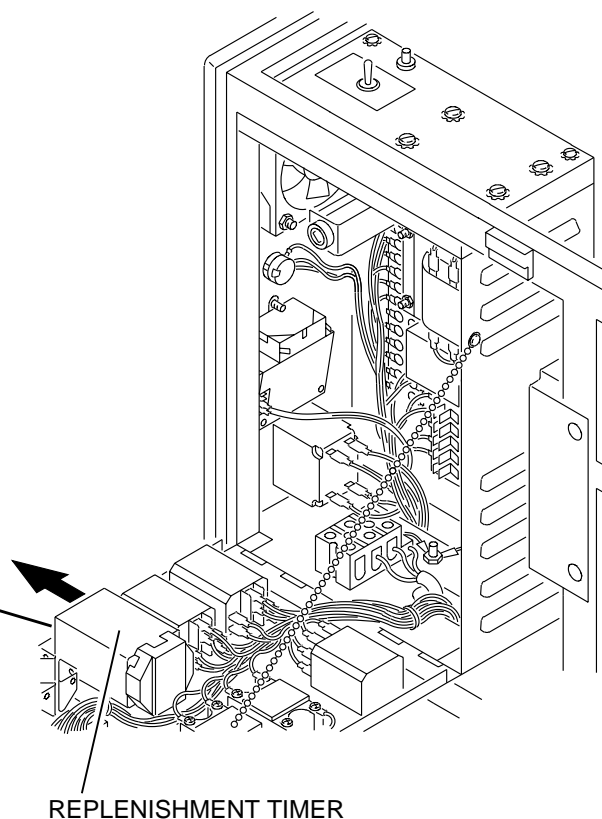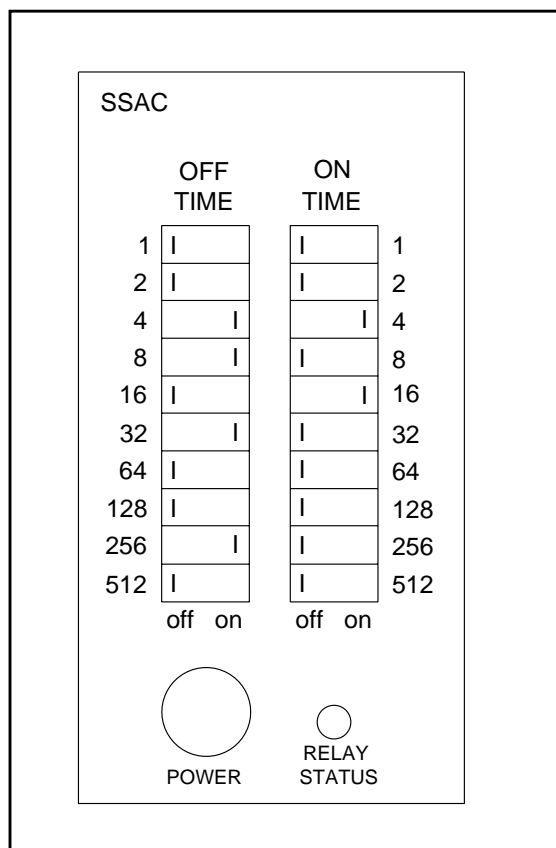

H112_0226HCA
H112_0226HC

- [1] Remove and discard the existing REPLENISHMENT TIMER.
- [2] Verify that the new REPLENISHMENT TIMER is set to:

"Off Time"	"On Time"
300 seconds	20 seconds

Note

The numbers on the REPLENISHMENT TIMER represent seconds.

- [3] Install the new REPLENISHMENT TIMER.

Publication History

Date	ECO No.	Pub. No.	Affected Pages	File Name	Notes
12FEB99	2590-355	3H9028	All	ii3408_1_12feb99.fm	First Printing
29APR99	2590-355	3H9028	All	ii3408_1_29apr99.fm	Revise

Kodak, Kodak Professional, Industrex, Industrial, and X-Omat are trademarks.

